



# **Foldlite One User Manual**

## Urban Fox Foldlite One

Overall Length	720 mm (28 in.)
Overall Width	640 mm (25 in.)
Folded Width	270 mm (11 in.)
Handle Height (Max)	880 mm (35 in.)
Handle Height (Min)	780 mm (31 in.)
Seat Width	480 mm (19 in.)
Seat Height	560 mm (22 in.)
Seat Depth	200 mm (8 in.)
User Weight (Max)	21 st. (135 kg)
Total Weight	5.45 kg
Wheels	200 mm (8 in.) Punchure Proof

### Parts Description




1. Hand Grip	2. Back Rest
3. Brake Lever	4. Pull Strap For Folding
5. Seat	6. Handle Tightening Knob
7. Rollator Bag	8. Rear Wheel
9. Stick holder	10. Front Wheel

## Assembly & User Instructions

### 1. Remove all packaging and set-out parts.


### 2. Inserting Front & Rear Wheels:

Install front and rear wheel sets as shown below.

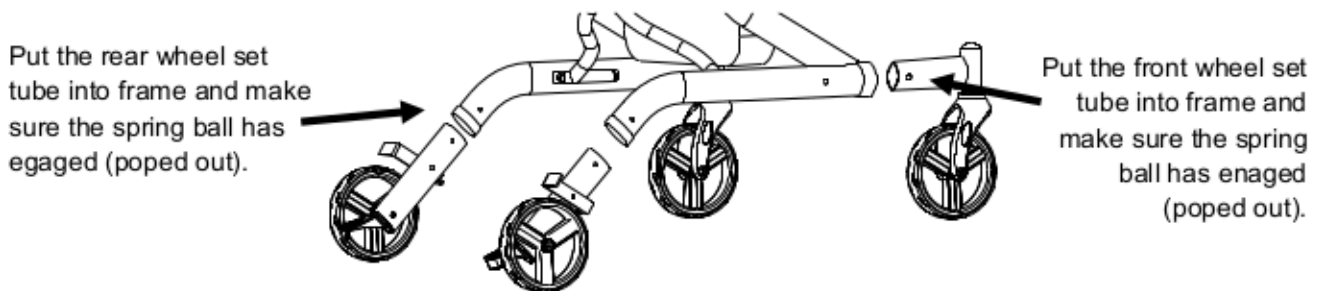

### 3. Opening Rollator:

To open the rollator, push down on either side of the cross bars as shown below. You should feel both sides 'click' into position.

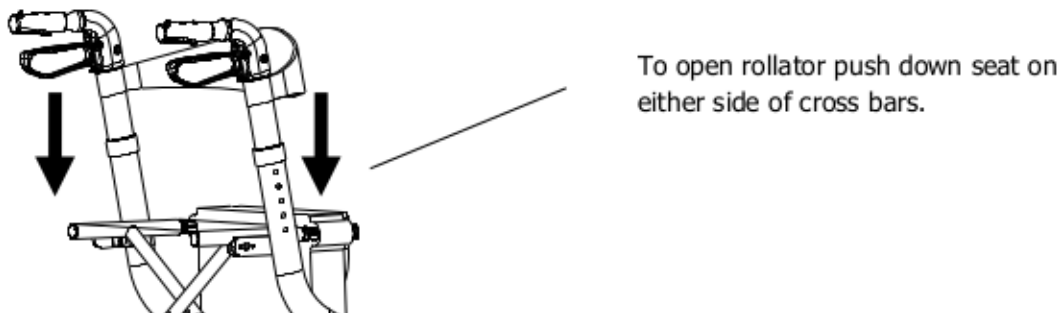


### **WARNING!**

**Be carefull not to trap fingers when pushing down seat.**


#### 4. Inserting the Handles:

Insert the Handles as shown and push down to the required height. Use the tightening knobs provided to secure in place at the desired height.

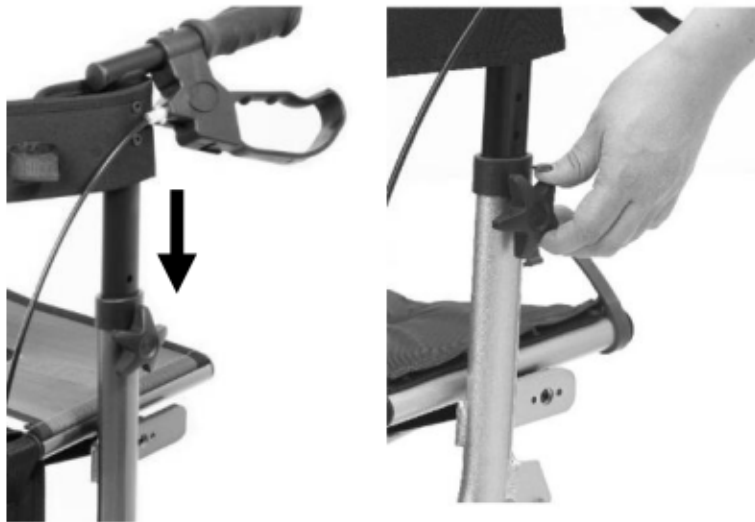

#### 5. Folding the Rollator

To fold the rollator pull up the strap in the middle of the seat as shown below.

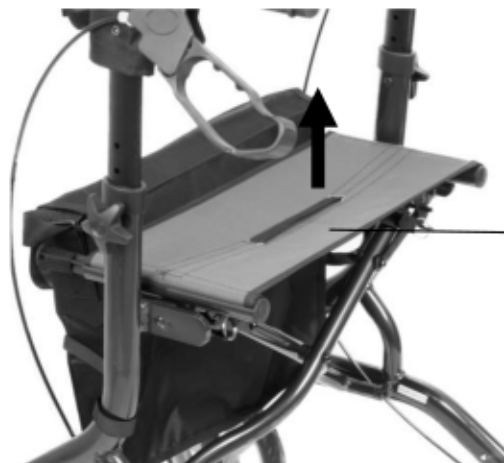


To fold rollator pull up the nylon strap in the center of the seat.

Note: As shown below a band is supplied to hold rollator together in folded position.



## 6. Rollator Bag:

Attach the bag to the rollator by placing the bag hoops over the two bag hangers positioned at the front of the seat as shown. Place the two velcro bag straps around the rollator frame to secure the lower section of the bag.

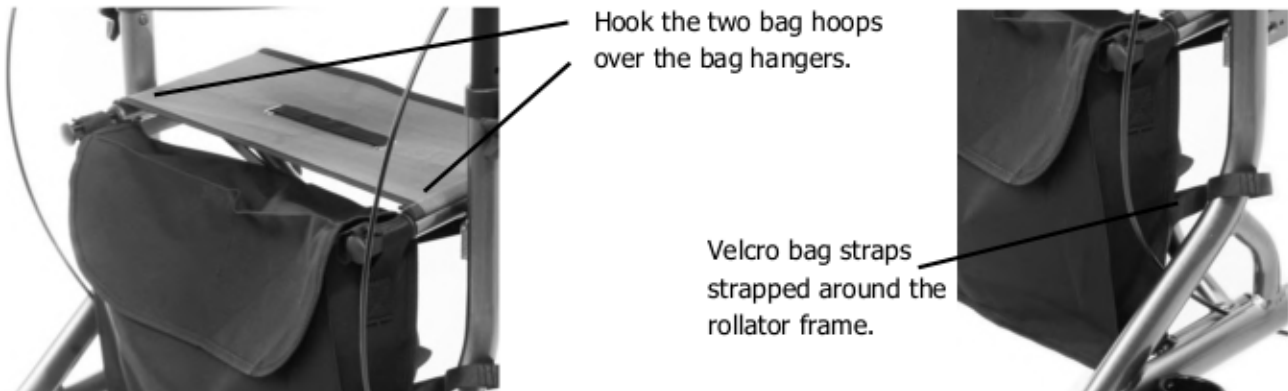

## 7. Back Strap:

Before using the rollator check that the back strap frame bolts are securely tightened and that the velcro straps are secure.



## 8. Stick Holder:

Attached to the rollator is a stick holder. Place your stick into the stick holder as shown and secure stick by using the velcro strap attached to the back rest.



### 9. Handle Height Adjustment:

To adjust the handle height up or down, undo the tightening knob and adjust the height. Tighten the knob back up making sure the handles are secure and level.



### 10. Braking:

For braking while moving, pull both brake levers upwards as shown and always brake slowly and steadily.



### 11. Locking the Brakes (Parking Position):

Press both brake levers down to apply and lock the brakes while the rollator is stationary.



**WARNING!**

**Regularly Test The Brakes**

## 12. Sitting on Rollator:

Before sitting lock the brakes into the parking position and make sure the rollator is on an even and firm surface and will not slide or topple over. Sit straight up and never lean backwards.



**WARNING! NEVER MOVE ROLLATOR WHILE SEATED**



**WARNING! ALWAYS APPLY PARKING BRAKES BEFORE SITTING**

### Maintenance and other

1. Do not use strong corrosive cleaners to clean or polish surface. This could damage rollator.
2. Store in dry location.
3. Be careful when folding and unfolding rollator to avoid injury by mechanism.
4. Always apply brakes when in stationary position or when seated on rollator.
5. Do not hang bags or other items from handles – this will add stress to the rollator and could also make it un-stable.

